



GENOSS

Peripheral DCB PTA Catheter

- Crystalline paclitaxel of $3.5 \mu\text{g}/\text{mm}^2$ coating with shellac and vitamin E

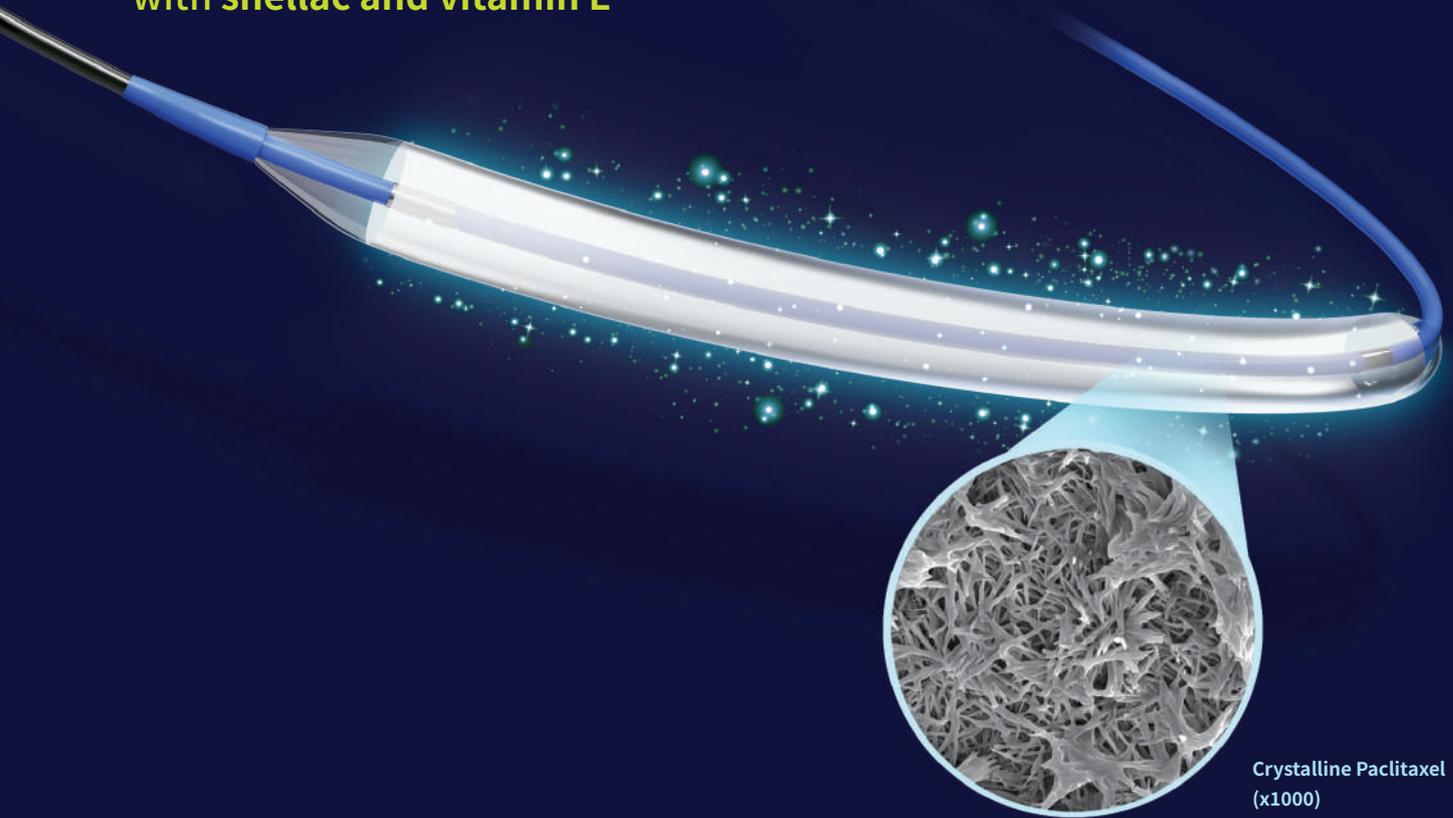


- Rapid deflation time by large (inflation & deflation) lumen diameter

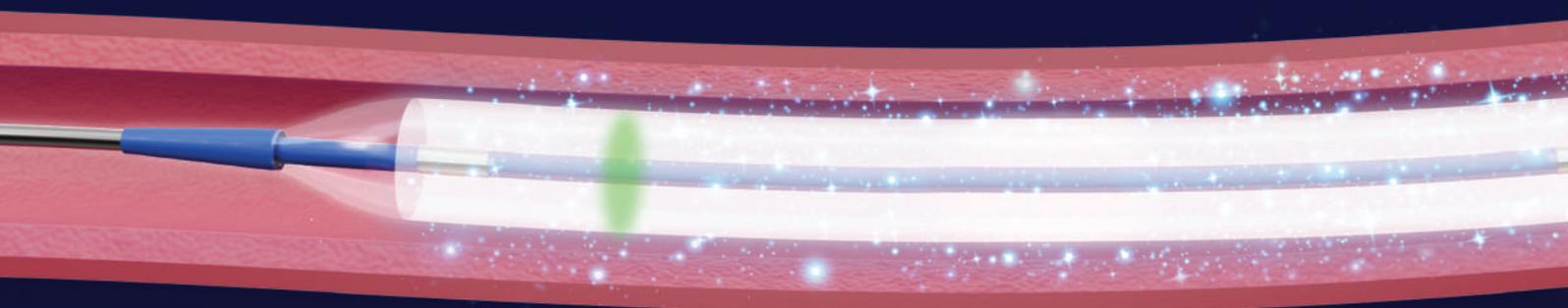


GENOSS Peripheral DCB PTA Catheter

Low drug loss by folding shape and covered surface rate

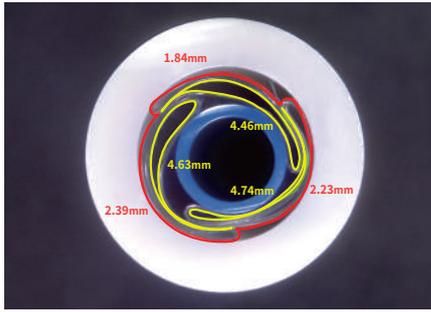

GENOSS

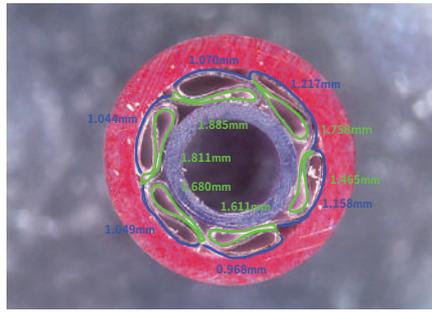

Company A


Sortation	Covered Balloon surface	Non-Covered Balloon surface	Covered surface rate
GENOSS (3-wing)	13.83 mm	6.46 mm	68.16 %
Company A (6-wing)	10.21 mm	6.50 mm	61.08 %

Comparison of 2 Paclitaxel-Coated Balloons

[J Vasc Interv Radiol 2025;1-8] Inclusion criteria: Rutherford classification 2-5, Lesion length ≤ 300mm


Clinical Outcomes

■ GENOSS ■ Company A

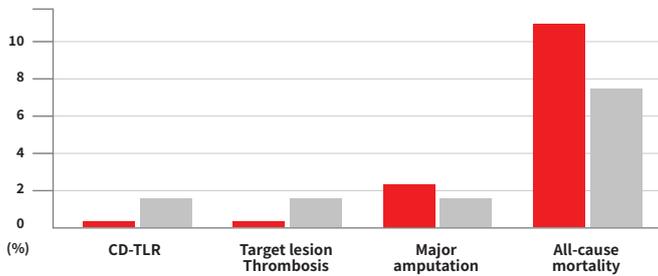

Angiographic Outcome

■ GENOSS ■ Company A

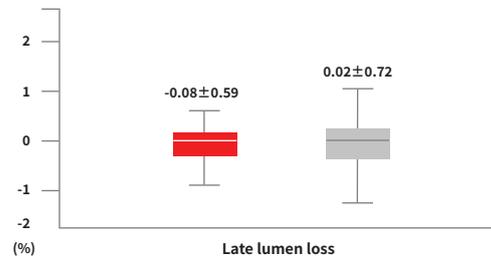

Order information

Applicable sheath: 6Fr (Catheter length : 50cm, 80cm,130cm)

Balloon Diameter	Balloon length	Product Codes	Balloon Diameter	Balloon length	Product Codes	Balloon Diameter	Balloon length	Product Codes	Balloon Diameter	Balloon length	Product Codes
Ø 4.0 mm	20 mm	PDEB-04-020-050	Ø 5.0 mm	20 mm	PDEB-05-020-050	Ø 6.0 mm	20 mm	PDEB-06-020-050	Ø 7.0 mm	20 mm	PDEB-07-020-050
		PDEB-04-020-080			PDEB-05-020-080			PDEB-06-020-080			PDEB-07-020-080
		PDEB-04-020-130			PDEB-05-020-130			PDEB-06-020-130			PDEB-07-020-130
	40 mm	PDEB-04-040-050		PDEB-05-040-050	PDEB-06-040-050		PDEB-07-040-050				
		PDEB-04-040-080		PDEB-05-040-080	PDEB-06-040-080		PDEB-07-040-080				
		PDEB-04-040-130		PDEB-05-040-130	PDEB-06-040-130		PDEB-07-040-130				
	60 mm	PDEB-04-060-050		PDEB-05-060-050	PDEB-06-060-050		PDEB-07-060-050				
		PDEB-04-060-080		PDEB-05-060-080	PDEB-06-060-080		PDEB-07-060-080				
		PDEB-04-060-130		PDEB-05-060-130	PDEB-06-060-130		PDEB-07-060-130				
	80 mm	PDEB-04-080-050		PDEB-05-080-050	PDEB-06-080-050		PDEB-07-080-050				
		PDEB-04-080-080		PDEB-05-080-080	PDEB-06-080-080		PDEB-07-080-080				
		PDEB-04-080-130		PDEB-05-080-130	PDEB-06-080-130		PDEB-07-080-130				
	100 mm	PDEB-04-100-050		PDEB-05-100-050	PDEB-06-100-050		PDEB-07-100-050				
		PDEB-04-100-080		PDEB-05-100-080	PDEB-06-100-080		PDEB-07-100-080				
		PDEB-04-100-130		PDEB-05-100-130	PDEB-06-100-130		PDEB-07-100-130				
	120 mm	PDEB-04-120-050		PDEB-05-120-050	PDEB-06-120-050		PDEB-07-120-050				
		PDEB-04-120-080		PDEB-05-120-080	PDEB-06-120-080		PDEB-07-120-080				
		PDEB-04-120-130		PDEB-05-120-130	PDEB-06-120-130		PDEB-07-120-130				
	150 mm	PDEB-04-150-050		PDEB-05-150-050	PDEB-06-150-050		PDEB-07-150-050				
		PDEB-04-150-080		PDEB-05-150-080	PDEB-06-150-080		PDEB-07-150-080				
		PDEB-04-150-130		PDEB-05-150-130	PDEB-06-150-130		PDEB-07-150-130				

Compliance chart

(atm)	(MPa)	Balloon diameter (mm)			
10 (NP)	1.01	Ø 4.0	Ø 5.0	Ø 6.0	Ø 7.0
11	1.11	4.05	5.04	6.07	7.06
12	1.22	4.10	5.08	6.12	7.12
13	1.32	4.15	5.12	6.18	7.15
14 (RBP)	1.42	4.20	5.16	6.24	7.19